



Vibrating Massager with Belt



Model: AB-751

Please read these instructions before use and keep for future reference.

Thank you for your purchase. Please read this user manual completely before use. Please follow the correct methods of usage, cleaning, maintenance and notes on safety. Please keep this user manual for future reference.

When using an electrical appliance, basic precautions should always be followed, including the following:


1. An appliance should never be left unattended when plugged in. Unplug from outlet when not in use and before putting on or taking off parts.
2. Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
3. Do not sleep with the unit operating. Long-time use may cause physical discomfort.
4. Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
5. Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
6. Never operate this appliance if it has a damaged cord or plug. If it is not working properly, if it has been dropped or damaged, or dropped into water. Contact Sunpentown for examination and repair.
7. Do not carry or drag this appliance by the power cord.
8. Do not sit or stand on the unit. Do not press hard against the unit.
9. Keep the cord away from heated surfaces.
10. Do not disassemble or transform the unit.
11. To disconnect, turn off unit and remove plug from the outlet.
12. Do not use within one hour before or after a meal.
13. Keep each session under 30 minutes.
14. After using this product, itchy skin is a normal phenomenon.
15. If anomalies or malfunctions occurs, please stop using immediately.
16. **THIS UNIT IS NOT A MEDICAL DEVICE.**

If you have the following conditions, please consult with your physician before using the unit.


1. Pregnant women or individuals who have undergone surgery and have inflamed or unhealed wounds.
2. Person who currently has an acute disease; has a history of heart problems; or has bone fractures that not healed completely; and other abnormal health conditions.
3. Person currently receiving treatment or with physical discomfort.
4. Cancer patients.
5. Spinal patients.
6. Children and those with intellectual disabilities must be under adult care when using.
7. While using this unit, if you start to feel dizzy, lightheaded or uncomfortable, please stop the unit at once and consult your physician. Do not continue use until you have the permission of your doctor.


- Do not use the massager in the bathroom or other damp/wet locations. Stop use immediately if unit gets wet.
- Please use the unit with an ambient temperature between 32°F-104°F
- Do not place unit near heaters and avoid sun exposure.
- Do not use outdoors. This unit is for household use only.

SAVE THESE INSTRUCTIONS

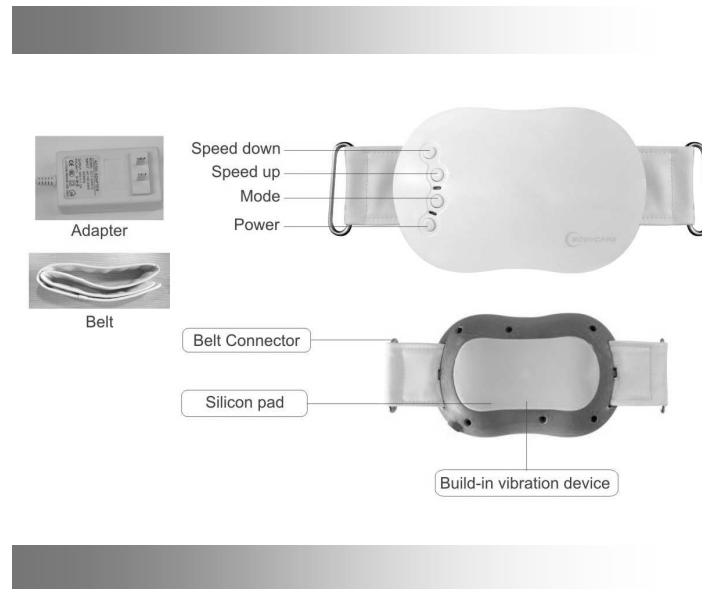

This massager can be used on various parts of the body, such as waist, hip, back and thigh.

1. Connect adapter to unit.
2. Plug adapter to power socket. A “beep” will be heard, unit is powered and in standby mode.
3. Strap belt through one end of the hoop. Place massager on desired body part, then strap belt around and through the other hoop. Secure in place.
4. Press the  button to begin. Unit will automatically operate under Intelligent mode.



5. Press  button to change to Manual mode.

6. To turn off, press  button again.

NOTE: Each session operates for 15 minutes and unit will shut-off automatically.

Mode Selection

The massager has two modes: Intelligent mode and Manual mode.

Default is Intelligent mode: power indicator and mode indicator are both lit. Press the  button and massager goes into manual mode (mode indicator goes off); pressing the  button again will return to intelligent mode and the cycle repeats.

Speed Option

Intelligent mode: speed cannot be changed. Unit cycles through different speeds and operates under a pre-programmed massage.

Manual mode: 5 intensity speeds. Press the **+** or **-** button to change the speed. Highest speed is 5 and lowest speed is 1. When switched over to manual mode, default is at speed 1.

NOTE: turn off unit and unplug from power when done.



Your Guarantee

If this product is found to be faulty as a result of faulty materials or workmanship within one year from date of purchase, it will be repaired free of charge.

This guarantee is subject to the following terms:

- Sunpentown must be notified of the fault.
- Proof of purchase must be presented to Sunpentown's nominated representative.
- The warranty will be void if the product is modified, misused or repaired by an unauthorized person.
- The warranty after repair will not be extended beyond the original one-year period.
- All replacement parts will be new or reconditioned.
- Parts, which are replaced, become the property of Sunpentown.
- The warranty applies for the use of the product in the USA only.

What is NOT COVERED:

- Warranty does not include freight charges.
- Incidental or consequential damage caused by possible defects with this product.
- Damage to product caused by improper power supply voltage, accident, fire, floods or acts of nature.
- Failure of product resulting from unauthorized modifications to the product.
- Improper installation or failure to perform the necessary maintenance.

This GUARANTEE is in addition to your Statutory Rights

SUNPENTOWN INTERNATIONAL INC.
14625 Clark Ave, City of Industry, CA 91745
Tel: 626-336-2299 · Fax: 626-336-7116
service@sunpentown.com
www.sunpentown.com

Cleaning

- Turn power off and unplug from power source before maintenance.
- Use a clean and dry cloth to wipe the unit.
- Never use thinner, benzene, or alcohol.

Storage

- For storage, it is suggested to place unit back in its original packaging and store in a cool, dry place.
- Do not disassemble.
- Keep out of the reach of small children.
- Do not store under direct sunlight or in places of high temperature.

Input	DC 12V / 1.8A
Rated power	15W
Universal adapter (UL)	100-240V / 50-60Hz
Massage unit dimension	7.13L x 4.63W x 2H in.
Belt	2.63W x 34L in.
Package dimension	12.25 x 8.25 x 3.63 in.
Net weight (massage unit only)	1.2 lbs
Gross weight	2.15 lbs