



## DEEP KNEADING FOOT MASSAGER

### OPERATING INSTRUCTIONS AB-763B




Thank you for your purchase. Please read this instructions manual thoroughly before use and keep in a safe place.


## FEATURES

- Elegant and ergonomic design.
- Infrared for heating treatment.
- Four groups of massager discs for comfortable kneading experience.
- Foot rollers with pressure nodes target vital reflex points on your soles, replicating the push and release technique used by professional reflexologists.
- Stimulates blood circulation, relieves tension and reduces stress.
- Detachable fabric covering for easy cleaning.

## NAME OF PARTS

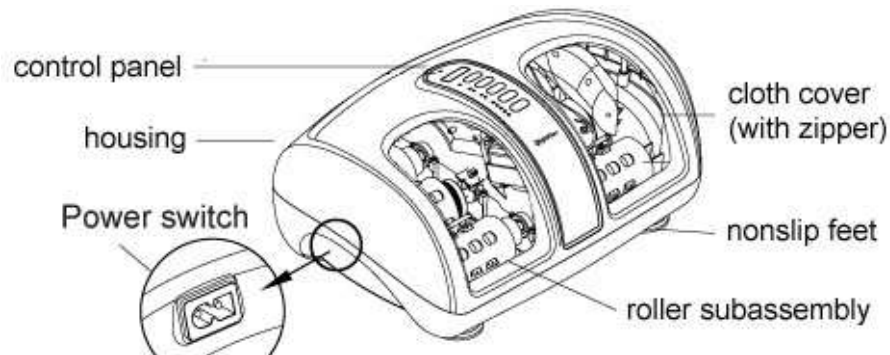

## CONTROL PANEL

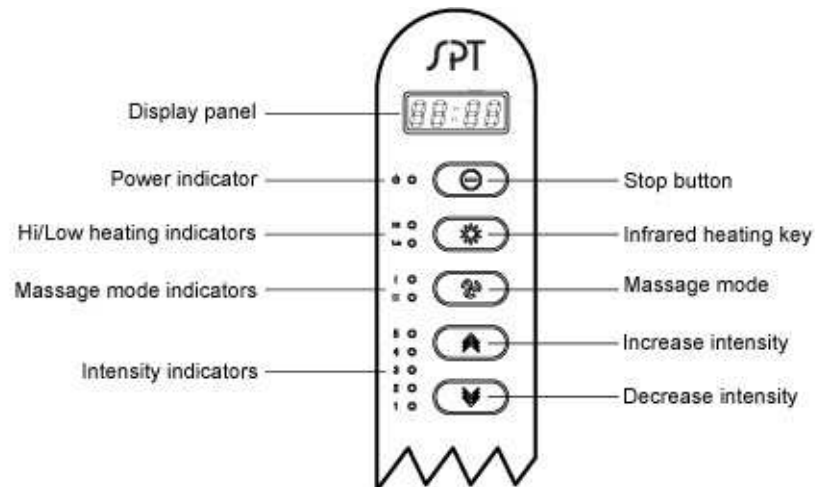

## DISPLAY PANEL

- Stand by mode: Power Indicator flashes and panel displays 00:00.
- When in session, panel alternately displays the remaining massage time and selected mode and intensity:  
12:58 = remaining session time  
HE-L = Heat Low / HE-H = Heat High  
P1 = Program 1 / P2 = Program 2  
.01 / .02 / .03 / .04 / .05 = corresponding intensity level

## OPERATING INSTRUCTIONS

- |  |  |  |   |
|--|--|--|---|
|  | Connect coupler terminal of power cord according to direction of socket. |  | Plug the unit into a properly grounded wall outlet.                           |
|  | Choose a chair of equal stature and place feet into the foot beds.       |  | Press to select High or Low infrared heating treatment.                       |
|  | Press to select message mode:<br>I = Program 1<br>II = Program 2         |  | Press to increase or decrease message intensity:<br>1 = Lowest<br>5 = Highest |
|  | Foot rollers with pressure nodes stimulates acupuncture points.          |  | Press STOP button to end session at any time. Unplug when not in use.         |

## NOTE

- A 15-minute massage session is recommended. The unit will automatically stop after 15 minutes.
- The massage may tickle, which is normal and a sign of blood circulation improvement.
- **IMPORTANT:** Do not sit or stand on the unit or exert excessive pressure on the foot beds. This may cause overheating and damage to the motor.
- The kneading discs will not stop instantaneously; they will complete one final cycle before stopping and returning to "home position".

## CAUTION

Consult your physician before use if you:

- Have a serious medical condition such as osteoporosis, skin diseases, heart disease, etc.
- Have an electronic medical device in your body, such as a pacemaker.
- May be pregnant.
- Are receiving medical care or are feeling unwell.
- Wish to use the unit on parts of the body that have previously been treated due to injury or sickness, or have eczema.
- Are experiencing pain on any part of your body.
- Have joint dysfunction such as rheumatism, hammer toe and gout.

## CAUTION

Instances when the massager should not be used. In the following events, please switch the power off immediately and disconnect power source:

- If water is accidentally spilled on the massager. This may result in electrical shock.
- If severe or abnormal pain is felt during massage. Stop operation immediately and consult your doctor.
- If you detect any malfunction or any other abnormal condition during operation.
- If there is a power failure, injury may occur when power is restored unexpectedly.
- If there is lightning.

## SAFETY

- Before connecting power, check that the voltage corresponds with the specification indicated for the massager.
- Do not disconnect the AC plug from the electrical outlet roughly.
- Do not connect or disconnect the AC plug from the electrical outlet with wet hands. May result in electrical shock or machine failure.
- When disconnecting from outlet, pull by grasping the plug and not the cord.
- Do not place the AC cord under the massager or other heavy objects.
- Do not wind the AC cord around the massager as it may damage the cord and result in fire or electrical shock.
- Do not operate the massager with a damaged cord or tamper with the cord. Contact Sunpentown immediately if cord is damaged.
- Do not use any attachments with this unit.
- Do not use the massager if the electrical outlet is excessively loose.
- Do not use in an excessively humid environment such as the bathroom.
- Do not use massager with wet feet. This may result in a short circuit or electric shock.
- Do not use the massager for any purpose, other than massage.

- Do not use the massager when the fabric cover is not attached or torn.
- Do not let children, or persons who are not capable of handling equipment, operate the massager.
- Please keep the massager out of reach of children.
- Do not insert metallic objects, or flammable materials into the main body. May result in electrocution, short circuit and possible fire.
- Do not use the massager within an hour after consuming a meal.
- Do not cover the massager during operation. May cause motor to overheat and subsequently breakdown.
- Do not unplug massager during operation, use the Stop button.

## PRODUCT CARE AND STORAGE

- Be sure to remove the AC plug from the electrical outlet before cleaning the massager.
- Do not use benzene, thinner or any other abrasive solvent on your massager.
- Use a dry cloth when cleaning the housing. Never use a damp cloth to clean the massager, power socket or control panel.
- When necessary, detach and wash the fabric cover for good hygiene.
- When not in use for a long period of time, it is advisable that the AC power cord be coiled and the massager stored in a dust and moisture free environment.
- Do not store the massager near intense heat or open flame. Do not leave the massager under direct sunlight for long period of time.

## ASSEMBLY AND REPAIR

The massager must only be serviced or repaired by Sunpentown Int'l Inc. Users must not disassemble or repair the massager themselves.

MODEL	AB-763B
DESCRIPTION	DEEP KNEADING FOOT MASSAGER
VOLTAGE	DC24V
POWER	60W
CYCLE	15 MINTUES



## Your Guarantee

If this product is found to be faulty as a result of faulty materials or workmanship within one year from date of purchase, it will be repaired free of charge.

This guarantee is subject to the following terms:

- Sunpentown must be notified of the fault.
- Proof of purchase must be presented to Sunpentown's nominated representative.
- The warranty will be void if the product is modified, misused or repaired by an unauthorized person.
- The warranty after repair will not be extended beyond the original one-year period.
- All replacement parts will be new or reconditioned.
- Parts, which are replaced, become the property of Sunpentown.
- The warranty applies for the use of the product in the USA only.

### What is NOT COVERED:

- Warranty does not include freight charges.
- Incidental or consequential damage caused by possible defects with this product.
- Damage to product caused by improper power supply voltage, accident, fire, floods or acts of nature.
- Failure of product resulting from unauthorized modifications to the product.
- Improper installation or failure to perform the necessary maintenance.

This GUARANTEE is in addition to your Statutory Rights

**SUNPENTOWN INTERNATIONAL INC.**  
**14625 Clark Ave. City of Industry, CA 91745**  
**Tel: 626-336-2299**  
**service@sunpentown.com**  
[www.sunpentown.com](http://www.sunpentown.com)